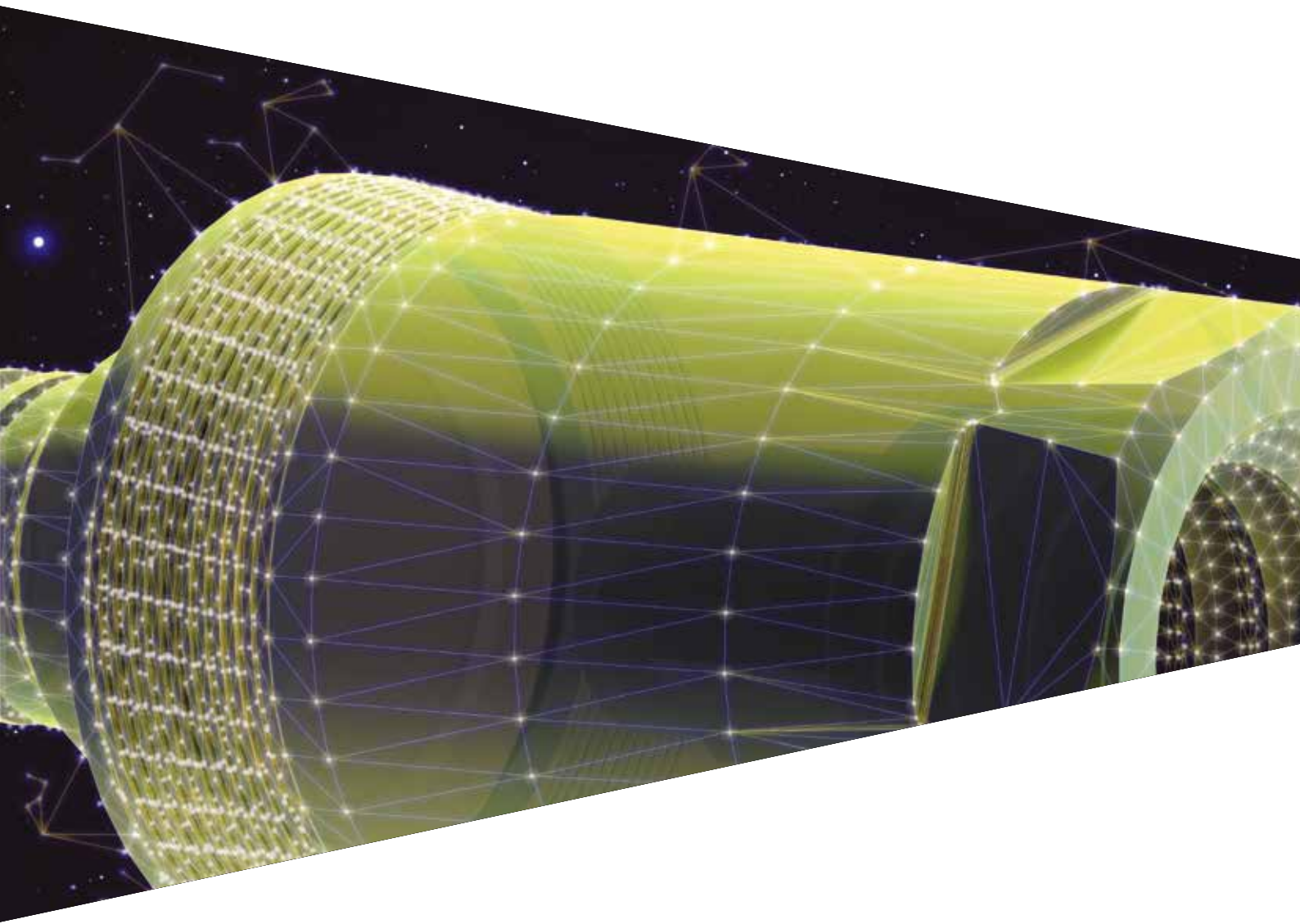


# XA Prosthesis

For SHELTA implants



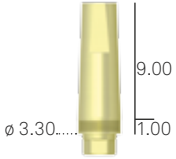
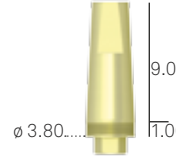
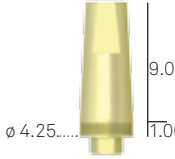
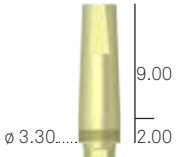
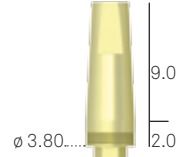
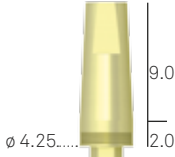






# XA posts for cemented prosthesis

These posts, made of Gr. 5 titanium, allow cemented prosthesis to be fitted using the One-Abutment-One-Time technique, in which the XA post is left screwed into the patient's mouth, reproducing its form and position in the laboratory model with great precision, using a transfer and a dedicated analog.

Posts are available with diameters of 3.30, 3.80 and 4.25 to permit Platform Switching on the three platforms of the Shelta implant.



## XA posts for cemented prosthesis

description	post $\varnothing$ 3.30	post $\varnothing$ 3.80 mm	post $\varnothing$ 4.25
for implants	Shelta 3.80	Shelta 3.80 - 4.25 - 5.00	Shelta 5.00
Premade XA post Repositionable Transgingival H. 1.00 mm	<b>SH-MD-F-330-1</b>  $\varnothing$ 3.30..... 9.00 1.00	<b>SH-MD-F-380-1</b>  $\varnothing$ 3.80..... 9.00 1.00	<b>SH-MD-F-425-1</b>  $\varnothing$ 4.25..... 9.00 1.00
Premade XA post Repositionable Transgingival H. 2.00 mm	<b>SH-MD-F-330-2</b>  $\varnothing$ 3.30..... 9.00 2.00	<b>SH-MD-F-380-2</b>  $\varnothing$ 3.80..... 9.00 2.00	<b>SH-MD-F-425-2</b>  $\varnothing$ 4.25..... 9.00 2.00
Single pack Pack of 10 pieces  Fixation screw with conical support	<b>L-VMS-180</b> <b>L-VMS-180-10</b>  M1.8.....	Use L-VMS-180	Use L-VMS-180
Analog for premade XA posts for cemented prosthesis	<b>SH-ANA-MD-F-330</b>  $\varnothing$ 3.30..... 9.00	<b>SH-ANA-MD-F-380</b>  $\varnothing$ 3.80..... 9.00	<b>SH-ANA-MD-F-425</b>  $\varnothing$ 4.25..... 9.00
Transfer for premade XA posts for cemented prosthesis	<b>SH-TRA-MD-F-330</b>  $\varnothing$ 3.30..... 10.00	<b>SH-TRA-MD-F-380</b>  $\varnothing$ 3.80..... 10.00	<b>SH-TRA-MD-F-425</b>  $\varnothing$ 4.25..... 10.00

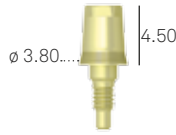
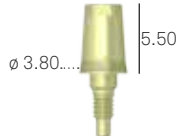
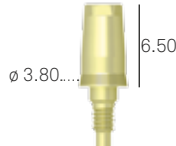

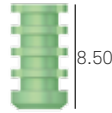


Recommended torque for definitive fastening of fixation screws: 20–25 Ncm.

# XA posts for screw-retained prosthesis

Intermediate XA abutments, with a feather-edge morphology, offer several prosthetic possibilities, from temporary through to final prostheses. As with the line of cemented prostheses, abutments for screw-retained prostheses can be fitted using the **One-Abutment-One-Time** technique, leaving the abutments screwed into the patient's mouth and reproducing their form and position in the laboratory model using a suitable analog. The only available diameter of 3.80 mm permits the greatest possible prosthetic simplification, and ensures Platform Switching.

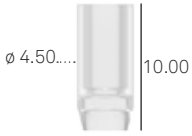
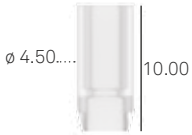
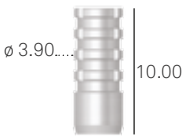
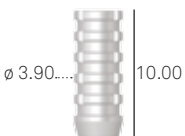



## XA posts for screw-retained prosthesis

∅ prosthetic component	∅ 3.80
for implants	Shelta 3.80 - 5.00
<p>Intermediate XA abutment H. 4.50 mm</p> <p>always combined with titanium sleeve when used for single-unit loading (SH-CTABU-F-380; SH-CTABU-F-380-ROT)</p>	<p><b>SH-ABU-F-TS-380-4</b></p> 
<p>Intermediate XA abutment H. 5.50 mm</p> <p>always combined with titanium sleeve when used for single-unit loading (SH-CTABU-F-380; SH-CTABU-F-380-ROT)</p>	<p><b>SH-ABU-F-TS-380-5</b></p> 
<p>Intermediate XA abutment H. 6.50 mm</p> <p>always combined with titanium sleeve when used for single-unit loading (SH-CTABU-F-380; SH-CTABU-F-380-ROT)</p>	<p><b>SH-ABU-F-TS-380-6</b></p> 
<p>Analog for intermediate XA abutments</p>	<p><b>SH-ANABU-F-380</b></p> 
<p>Transfer for intermediate XA abutments</p>	<p><b>SH-TRABU-F-380</b></p> 
<p>Transfer screw</p>	<p><b>SH-VTRABU-F-200</b></p> 
<p>Healing cap in PEEK for intermediate XA abutments</p>	<p><b>SH-CG-ABU-F-380</b></p> 

Recommended torque for definitive fastening of directly screwed posts: 25–30 Ncm.

## Accessories for screw-retained prosthesis on XA posts

for implants	Shelta 3.80 - 5.00
Castable sleeve for XA abutments Repositionable Fixation screw included	<b>SH-CCABU-F-380</b> 
Castable sleeve for XA abutments Non-repositionable Fixation screw included	<b>SH-CCABU-F-380-ROT</b> 
Titanium sleeve for XA abutments Repositionable Fixation screw included	<b>SH-CTABU-F-380</b> 
Titanium sleeve for XA abutments Non-repositionable Fixation screw included	<b>SH-CTABU-F-380-ROT</b> 
Fixation screw for XA abutments	<b>A-PLAIN-VP200</b> 

Recommended torque for definitive fastening of directly screwed structures on posts: 20–25 Ncm.





rev.07-16  
v. 01



**Sweden & Martina S.p.A.**

Via Veneto, 10  
35020 Due Carrare (PD), Italy  
Tel. +39.049.9124300  
Fax +39.049.9124290  
info@sweden-martina.com  
[www.sweden-martina.com](http://www.sweden-martina.com)

**Sweden & Martina Mediterranea S.L.**

Sorolla Center, Oficina 801  
Avda Cortes Valencianas 58, 8pl  
46015-Valencia, Spain  
Tel. +34.96.3525895  
Tel. 900535617  
info.es@sweden-martina.com

**Sweden & Martina Ltd**

Unit 45  
Basepoint Business Centre  
Metcalf Way, Crawley,  
West Sussex, RH11 7XX, UK  
Toll free 0800 1123575  
info.uk@sweden-martina.com

**Sweden & Martina Inc.**

c/o DCI Management  
301 Pleasant St.  
Abbottstown, 17301 PA, US  
Toll free 844-8MARTINA  
844-862-7846  
info.us@sweden-martina.com

The implants, prosthetic components and surgical instruments illustrated in this depliant are medical devices manufactured by Sweden & Martina SpA. The articles illustrated are compliant with the ISO 9001 and ISO 13485 standards, and are registered as CE Mark (Class I) and CE Mark 0476 (Class IIA and Class IIB) in accordance with the European Medical Devices Directive 93/42/EEC and European Directive 2007/47/EC.

Sweden & Martina production facilities manufacture medical devices in accordance with the cGMPs applicable in the USA and other countries.



Some products may not be available on all markets.

All brand names present in this depliant are the property of Sweden & Martina, except for those products for which other indications are given.

These products are intended for dental clinics and dental technology laboratories, and are not intended for direct sale to patients.

It is prohibited to sell, duplicate or disclose the products illustrated in this depliant without the prior written consent of Sweden & Martina S.p.A.

For further information on the products illustrated, including indications, contraindications, warnings, precautions and potential side-effects, please consult the Sweden & Martina S.p.A. website.

The contents of this depliant are correct at the time of publication. Please contact Sweden & Martina S.p.A. for information on any subsequent updates.